

# ELECTRODE PLACEMENT GUIDE

All about placement of  
Gel Electrode pads



Home BCR System  
SIGNATURE V2.0

## NOTE:

BCR-Therapie system used in Power2Cell device is a Holistic and Alternative therapy system and do not serve as a substitute for professional medical advice, examination, diagnosis or treatment.



### FACIAL BEAUTY

- A1 BEAUTY ACNE/SCARS
- A2 BEAUTY CELLULITE
- A3 FACE TONING
- A4 SKIN CONDITIONING
- A5 BEAUTY WRINKLES

### GENERAL HEALTH

- B1 ASTHMA
- B2 INFLAMMATION
- B3 SCARS
- B4 SINUS INFLAMMATION
- B5 SWELLING
- B6 TENSION
- B7 WOUNDS

### INFECTION

- C1 INFECTION - COLD
- C2 INFECTION - VIRUS

### INJURY HEALING

- D1 BRUISE
- D2 DISTORTION
- D3 FRACTURE
- D4 LIGAMENT INJURY
- D5 PULLED MUSCLE
- D6 TENDON IRRITATION
- D7 TENDON INJURY

### PAIN THERAPY

- E1 NECK PAIN SPINE
- E2 ARTHRITIS
- E3 CARPAL TUNNEL
- E4 FROZEN SHOULDER
- E5 HEEL PAIN
- E6 JOINT PAIN ACUTE
- E7 JOINT PAIN CHRONIC
- E8 MUSCLE PAIN
- E9 BACK PAIN SPINE

### PALSY REHAB

- F1 PALSY REHAB MUSCLE
- F2 PALSY REHAB NERVE
- F3 PALSY REHAB SPASTICITY
- F4 PALSY REHAB STIMULATION
- F5 ENERGY MODE

### SPORTS THERAPY

- G1 CONTUSION
- G2 MUSCLE RELAXANT
- G3 MUSCLE TREATMENT
- G4 TENNIS ELBOW
- G5 MUSCLE WARM UP
- G6 MUSCLE WARM DOWN

### WELLNESS

- H1 ALKALISATION
- H2 BLOOD CIRCULATION
- H3 DETOXIFICATION
- H4 LYMPH - FAT BURN
- H5 IMMUNE BOOST
- H6 LYMPH DRAINAGE
- H7 REJUVENATION
- H8 VEGETATIVE BALANCE

## Programs :

B2 INFLAMMATION for Headache

B6 TENSION for Headache & Migraine


F2 PALSY REHAB NERVE for Migraine, Parkinson's.

H8 VEGETATIVE BALANCE for Stress Relief


Indication : MIGRAINE / PARKINSON'S / STRESS RELIEF

### Option 1

Repeat reversing A and B

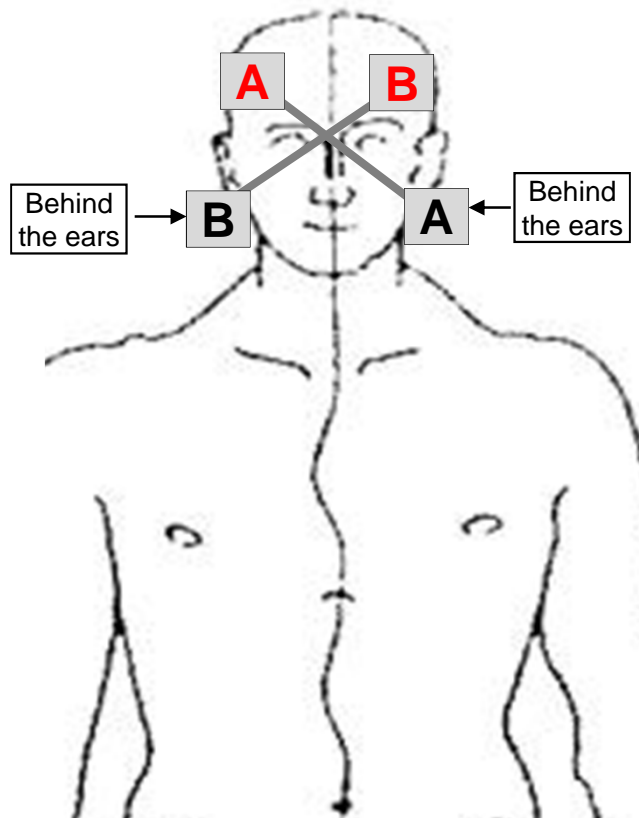

### Option 2

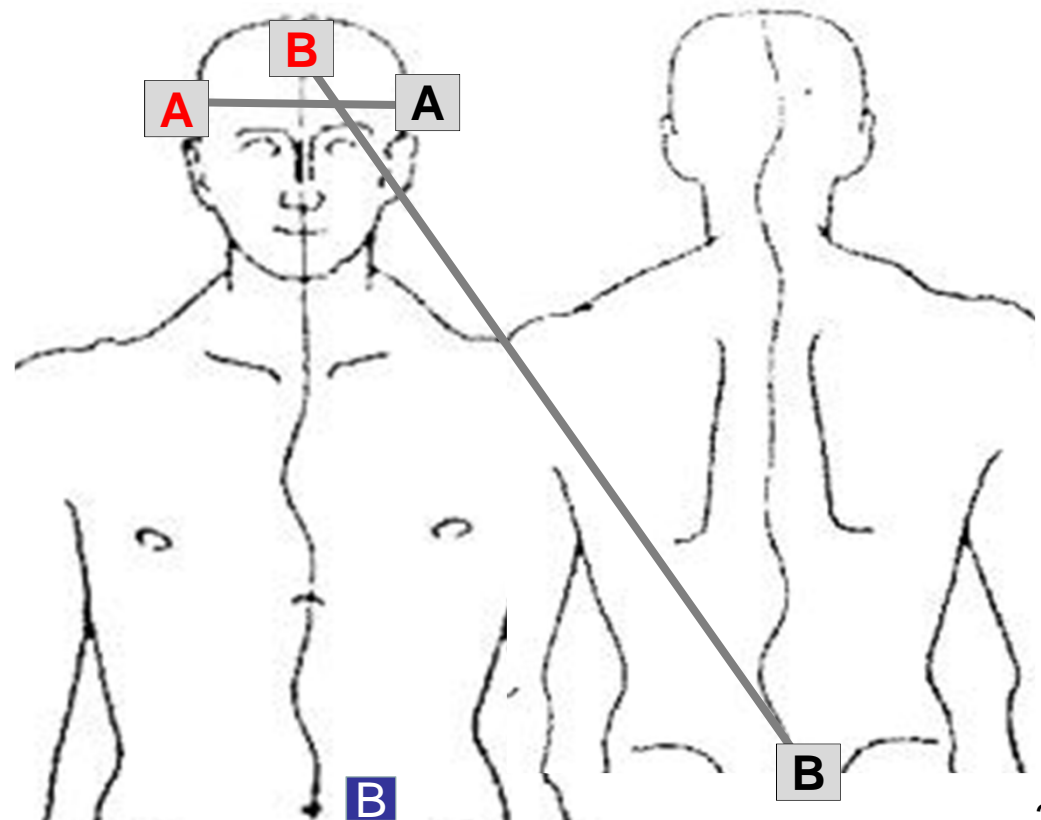


Indication : **NECK PAIN / OFFICE SYNDROME**

Programs :



**E1 NECK PAIN SPINE**

**E8 MUSCLE PAIN**

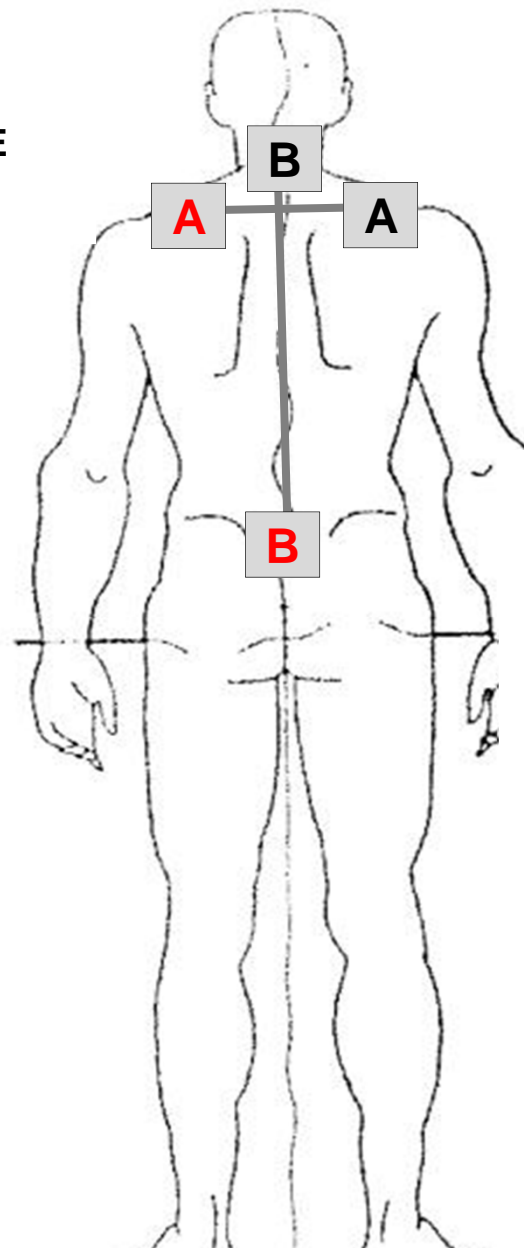

Place channel A  
horizontally across the  
target area


Include the arms as well (below)  
if the pain radiates to the arms

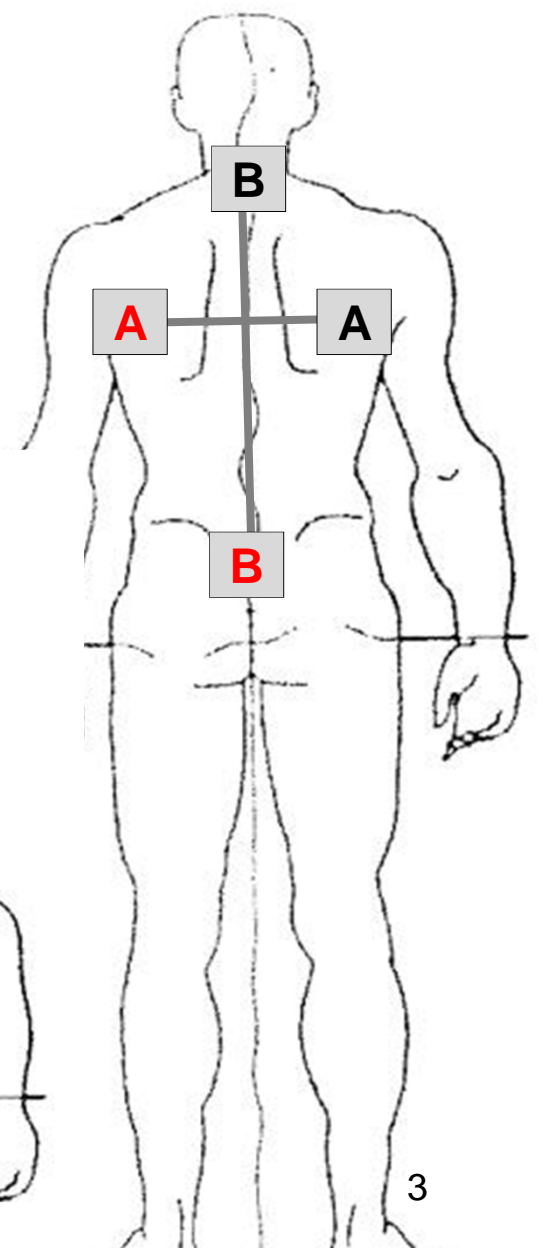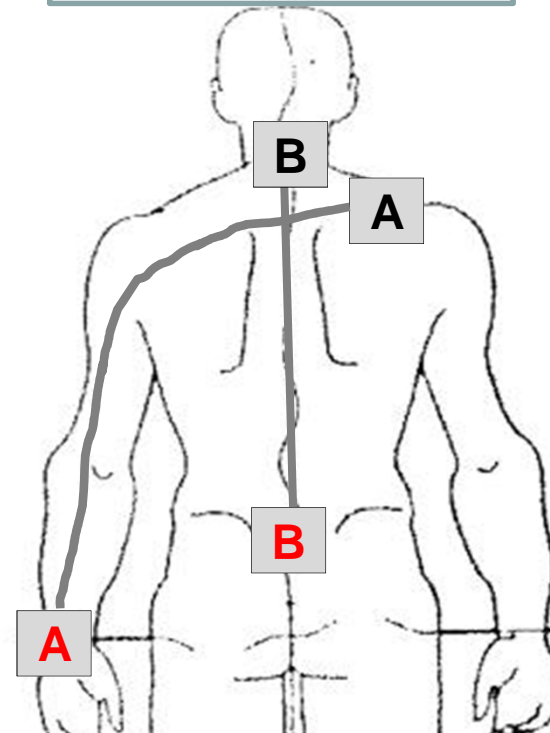

Indication : **FROZEN SHOULDER / SHOULDER JOINT PAIN**


Program :

**E4 FROZEN SHOULDER**

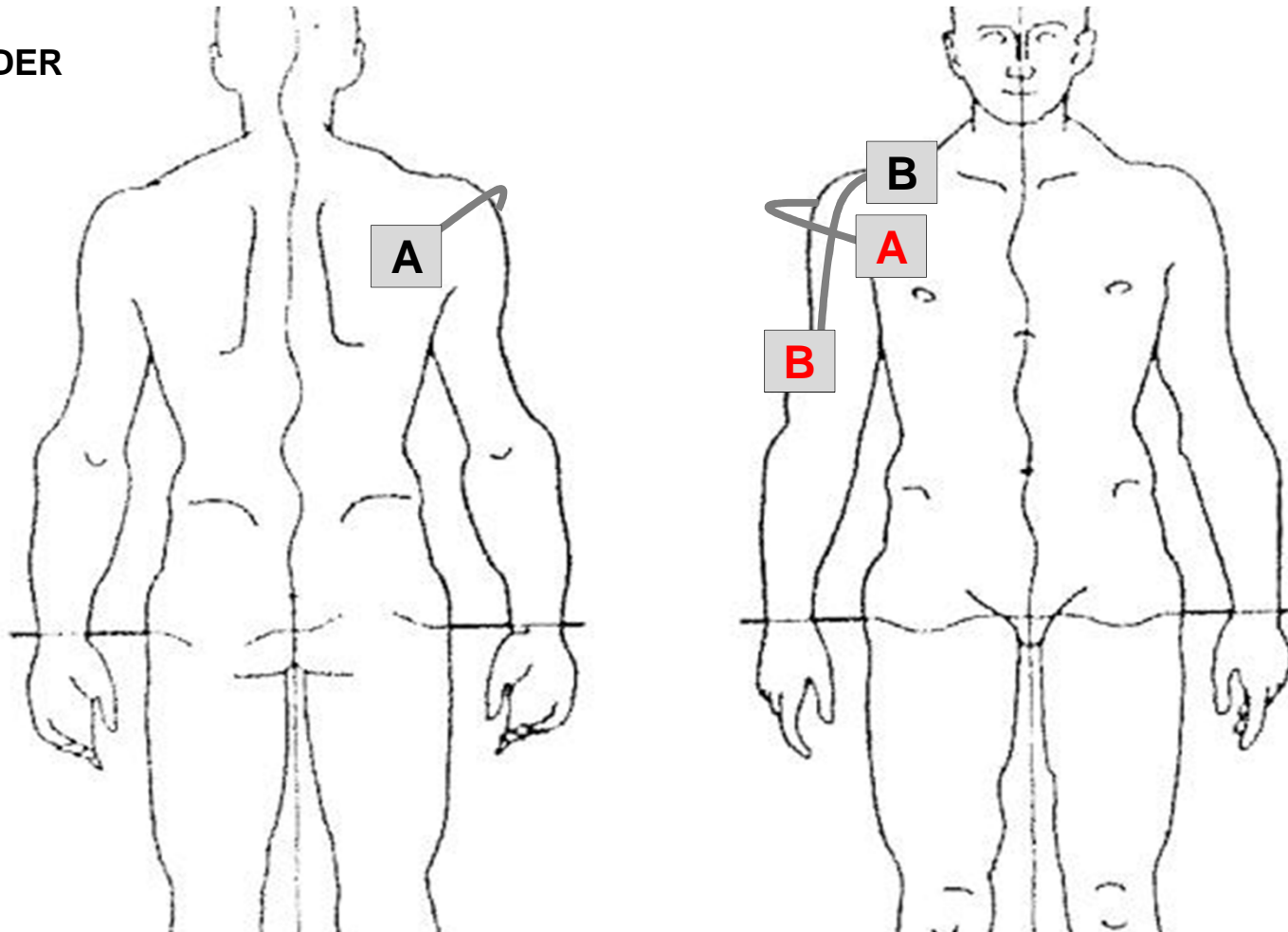

Indication : **TENNIS ELBOW / ELBOW JOINT PAIN**

Program : **G4 TENNIS ELBOW**

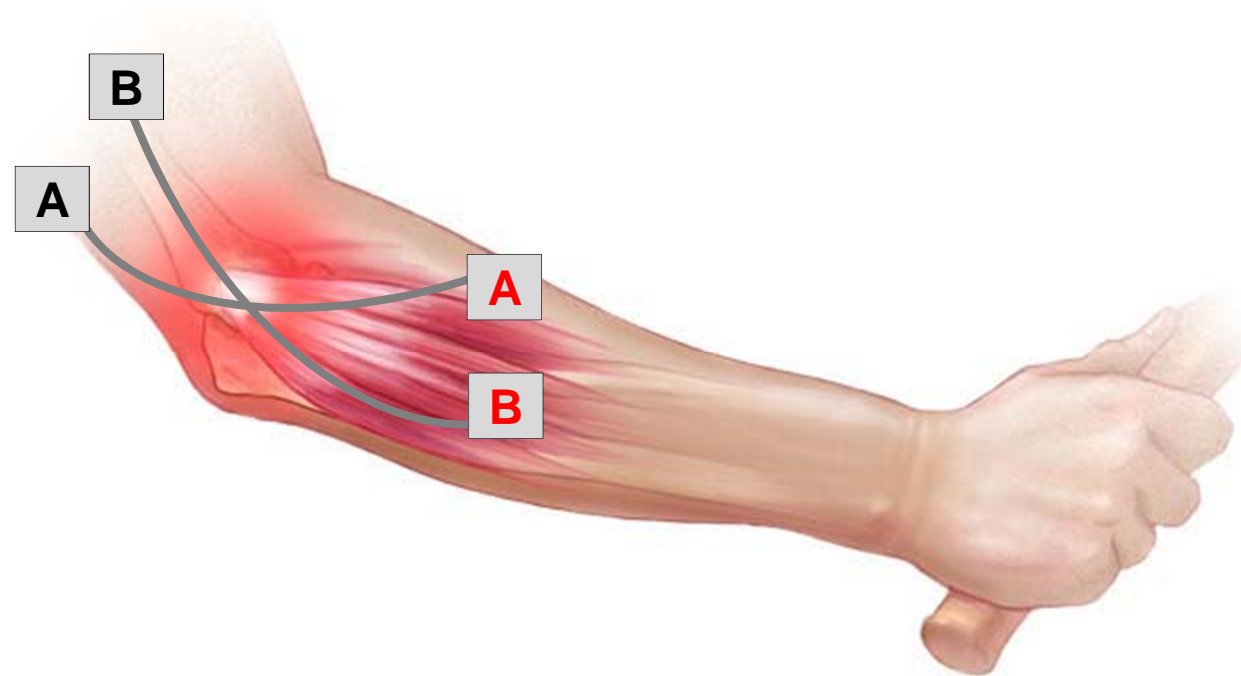


Indication : **BACK PAIN, SLIPPED DISC, SCIATICA**

Program :

**E9 BACK PAIN SPINE**

**G3 MUSCLE PAIN**

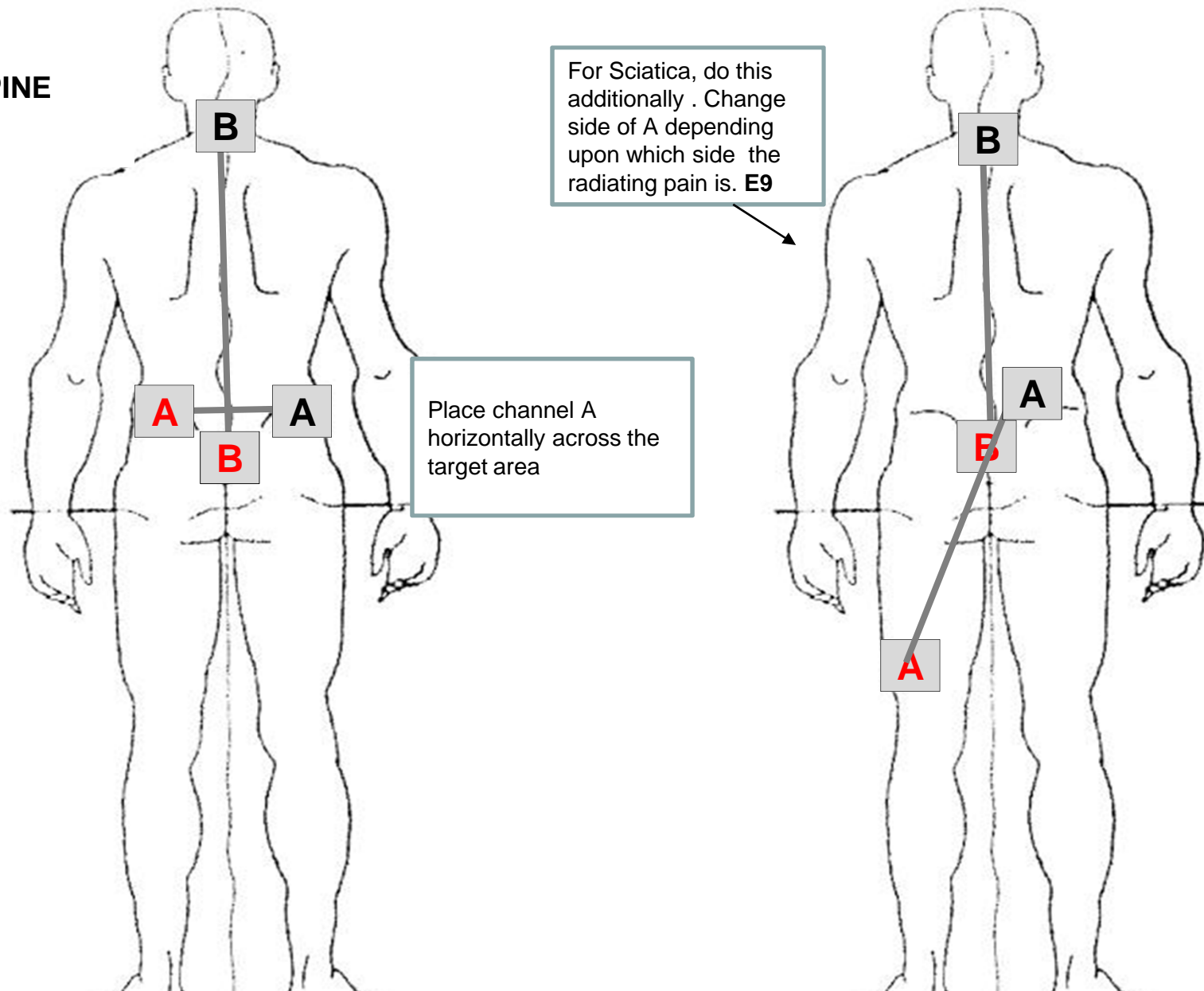


**Indication : BUTTOCK PAIN, HIP JOINT PAIN, AND SCIATICA CAUSED BY PIRIFORMIS MUSCLE**

**Program :**

**E9 BACK PAIN SPINE**

**G3 MUSCLE TREATMENT**

**E2 ARTHRITIS**

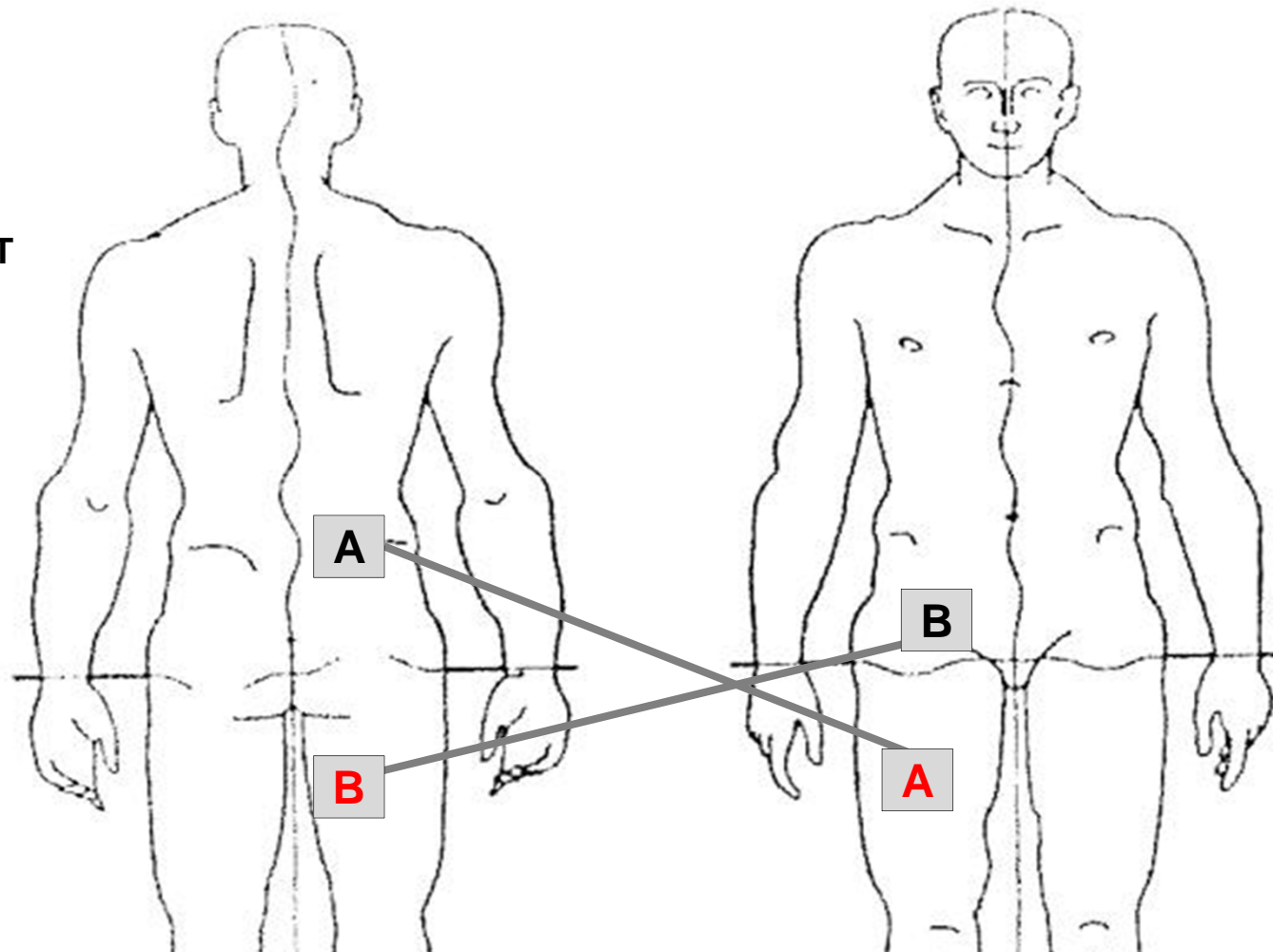

Indication : **KNEE PAIN**

Program :


E2 ARTHRITIS

E6 JOINT PAIN ACUTE


E7 JOINT PAIN CHRONIC

D4 LIGAMENT INJURY

Option 1

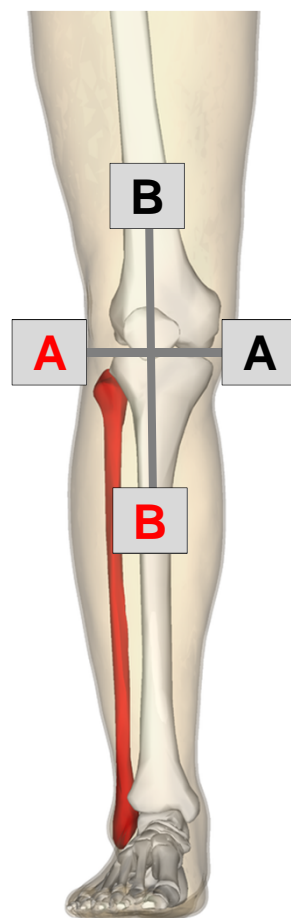

Option 2

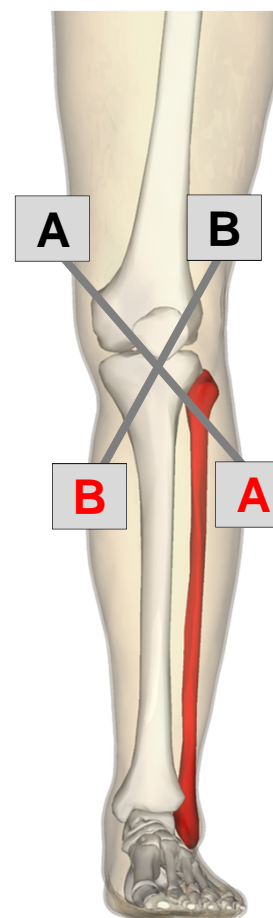




Indication : **SPRAINED ANKLE, ANKLE PAIN**

Program :

D2 DISTORTION

D6 TENDON IRRITATION

D7 TENDON INJURY

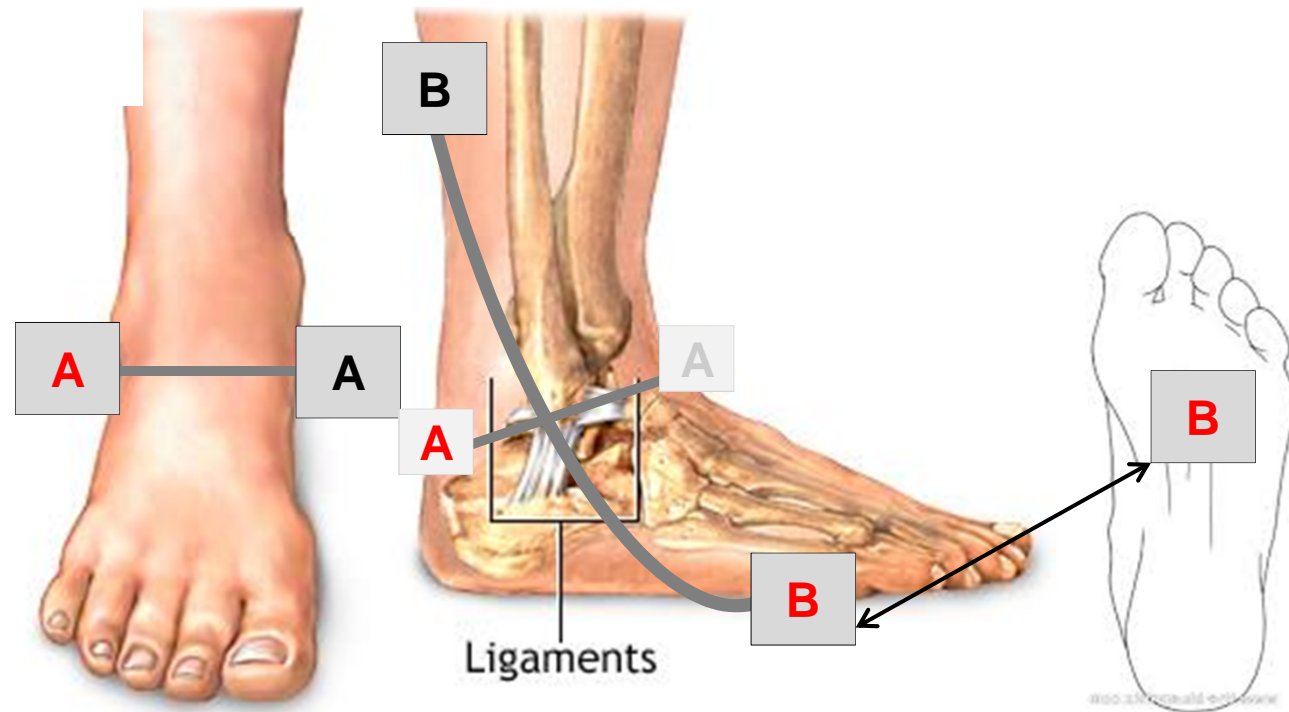


Indication : **ASTHMA , LUNG INFECTION**

Program :

**B1 ASTHMA**

**B2 INFLAMMATION**

**C1 INFECTION COLD**

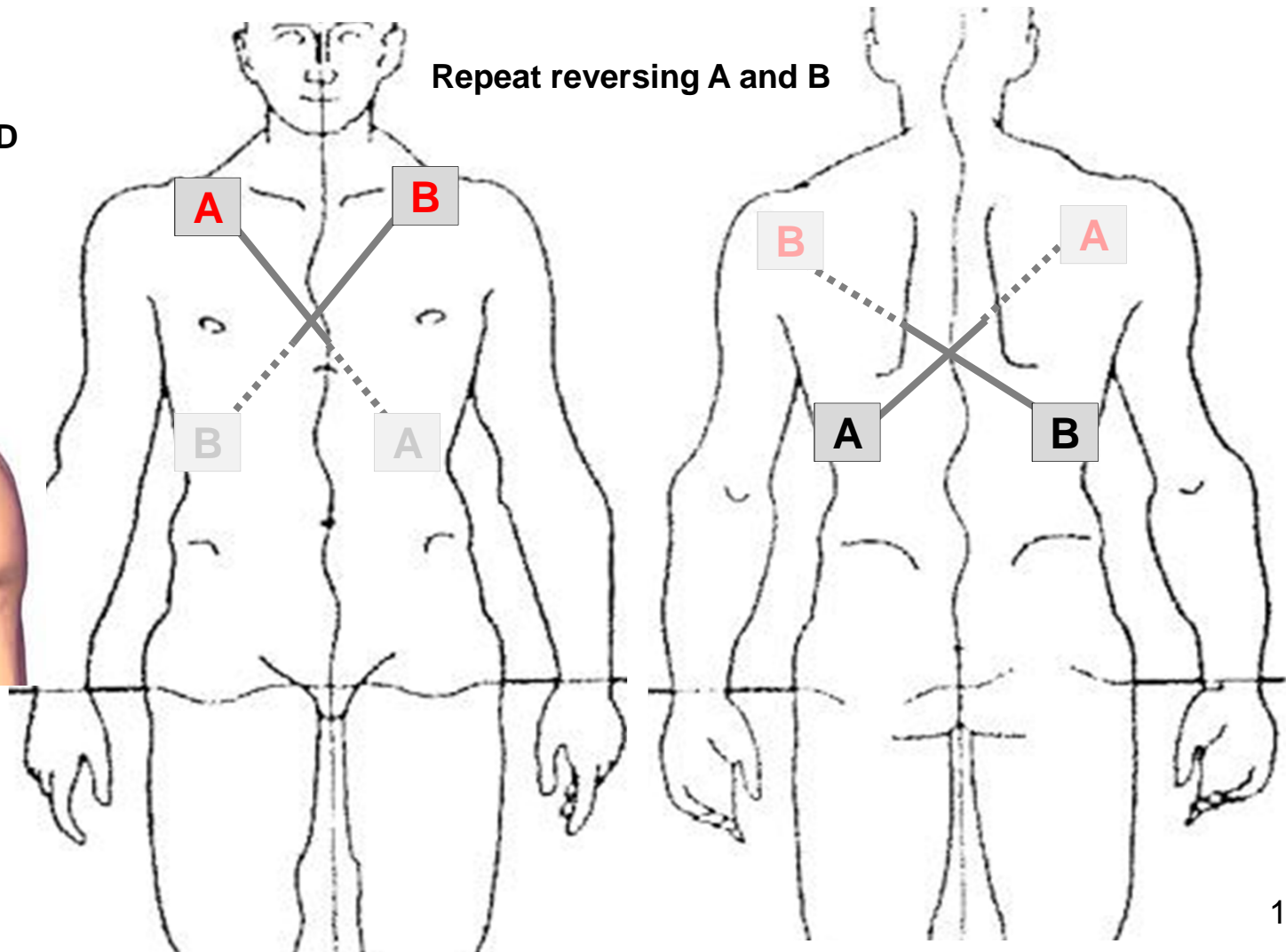

Indication : **MUSCLE WEAKNESS**


Program :

**H2 BLOOD CIRCULATION**

Repeat reversing A and B

**F4 PALSY REHAB STIMULATION**

Recommended for Erectile Dysfunction  
and Urinary Incontinence based on  
anecdotal evidences only, not  
supported by clinical studies yet.

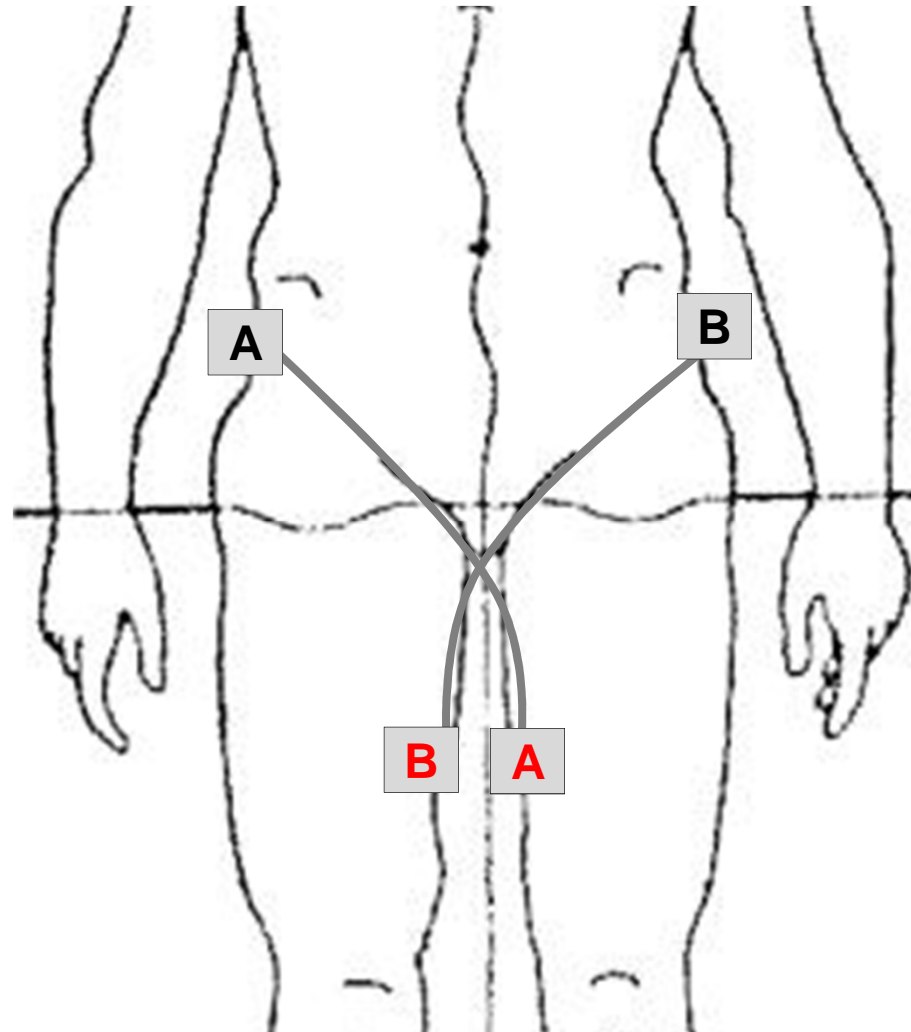


Indication : **SINUS INFLAMMATION**

Program :

**B4 SINUS INFLAMMATION**

Repeat reversing A and B

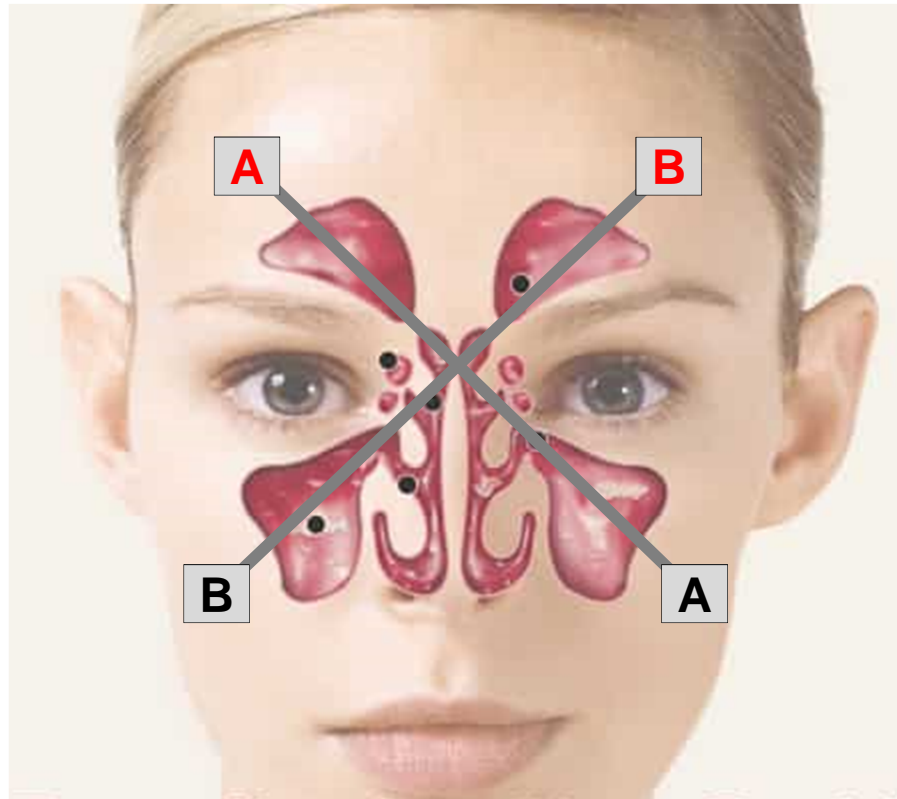


Indication : **WHOLE BODY (SYSTEMIC)**


Programs :

This is the systemic application. Can be used with all programs for systemic applications.

Step 1  
Repeat reversing  
A and B

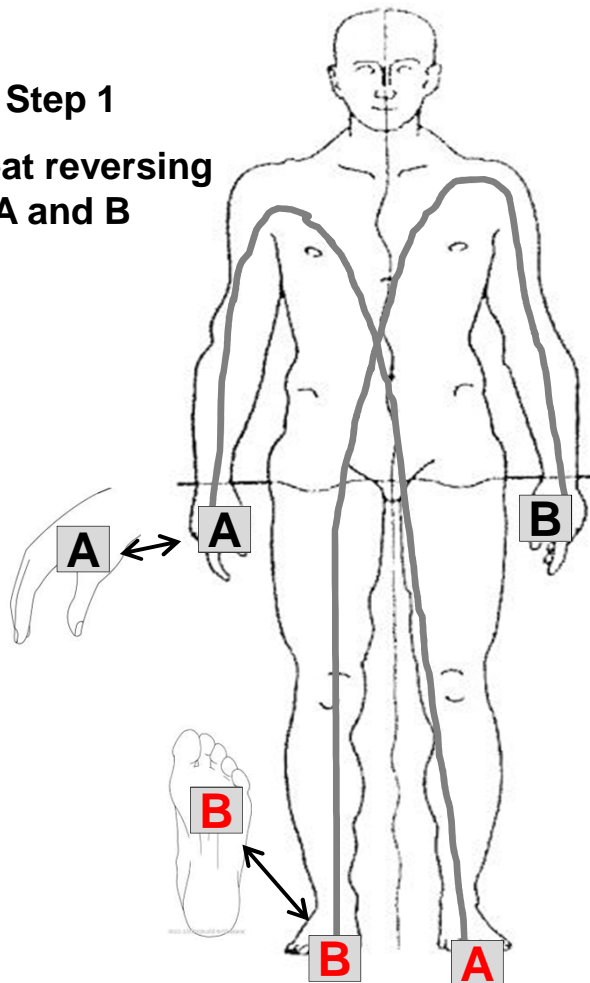

Step 2

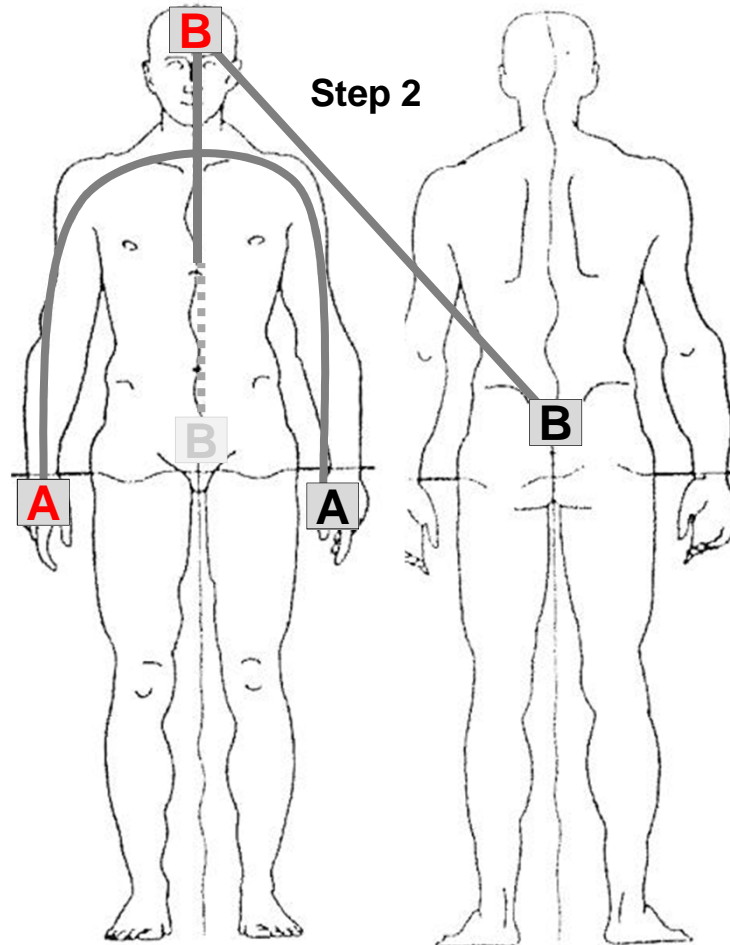


# STROKE RECOVERY / MUSCLE DYSTROPHY / SEMI-PARALYSIS

## Programs : F2 PALSY REHAB NERVE


Step 1 of 6

Repeat reversing A and B

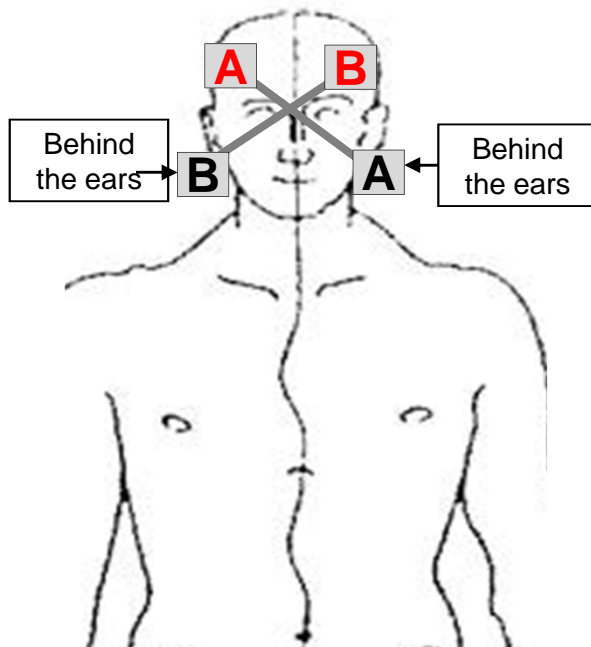

## Programs : F2 PALSY REHAB NERVE


Step 2 of 6

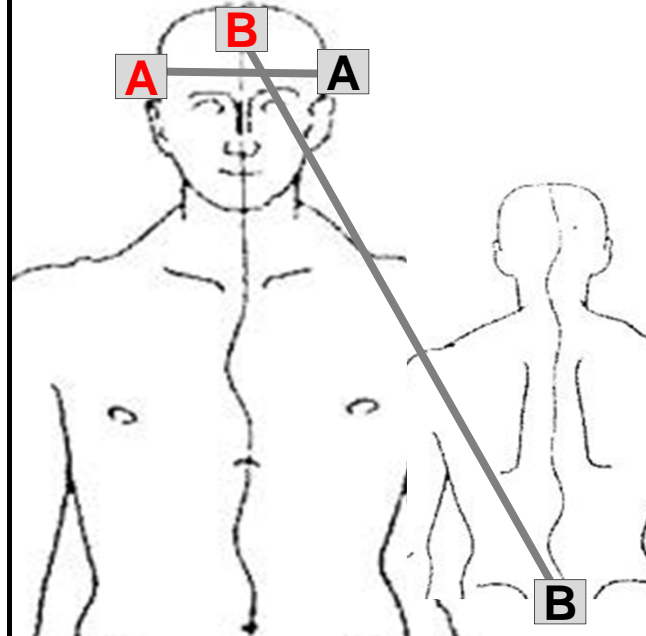

## Programs : F4 PALSY REHAB STIMULATION

Step 3 of 6


Repeat reversing A and B

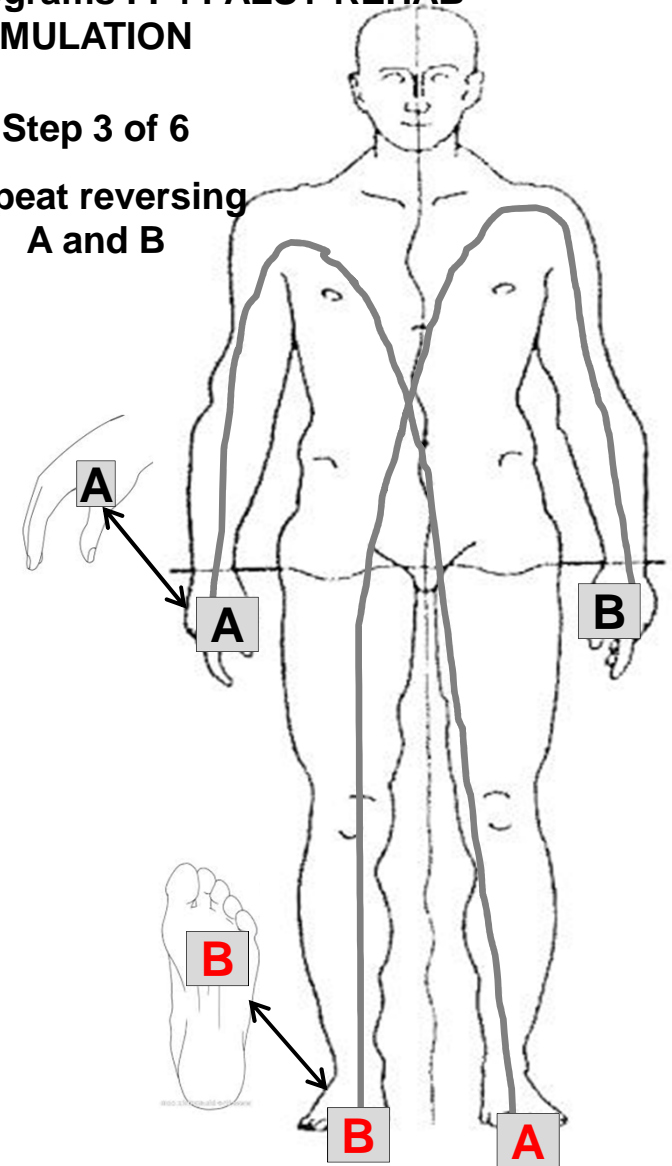



# STROKE RECOVERY / MUSCLE DYSTROPHY / SEMI-PARALYSIS

Programs : F2 PALSY REHAB NERVE, F4 PALSY REHAB STIMULATION

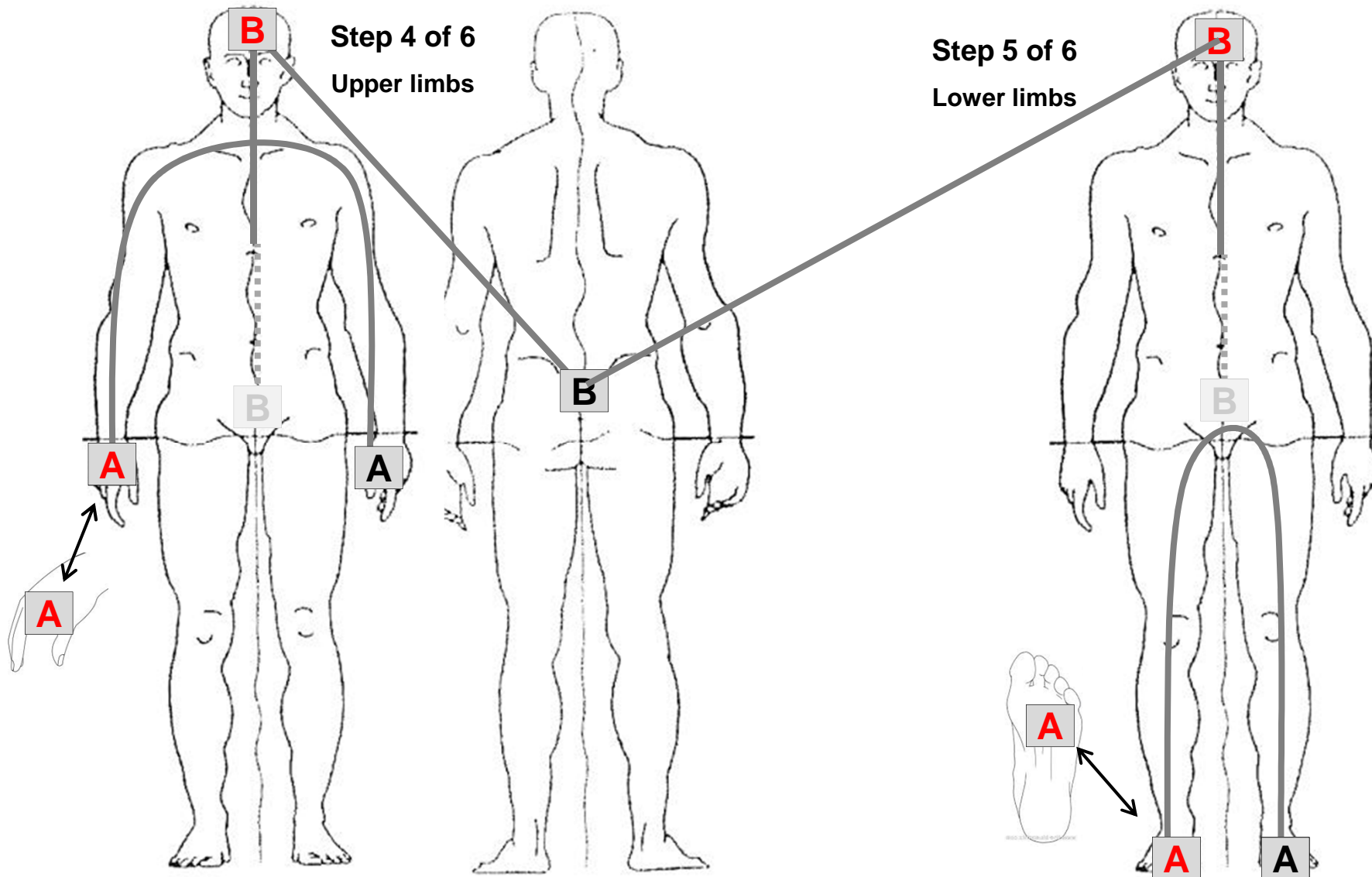

## STROKE RECOVERY / MUSCLE DYSTROPHY / SEMI-PARALYSIS

Step 6 of 6

**SPEECH / SWALLOWING**

Repeat reversing A and B

**Program:**

**F4 PALSY REHAB STIMULATION**

